EN-US Reprocessing instructions PL Instrukcja preparacji HU Előkészítési útmutató CS  
 Návod ke zpracování SK Návod na úpravu RO Instrucțiuni de pregătire HR Upute za  
 pripremu BG Инструкции за повторна обработка SL Navodila za pripravo SR Uputstvo za  
 održavanje LV Apstrādes instrukcija ET Töötlemisjuhend LT Apdorojimo instrukcija RU  
 Инструкции по обработке MK Упатство за повторна обработка TH คำแนะนำการประมวลผล  
 ใหม่ MS Arahan penyediaan VI Hướng dẫn xử lý AR لدلجة المع

LMT 67506h 01/2025 EN-US, PL, HU, CS, SK, RO, HR, BG, SL, SR,  
 LV, ET, LT, RU, MK, TH, MS, VI, AR

Nasal	Full Face	
JULIA	LENA	JOYCEone Full Face NV
CARA	LENA NV	JOYCEeasy Full Face
JOYCEone	CARA Full Face	JOYCEeasy next Full Face
JOYCEeasy	JOYCEone Full Face	JOYCEeasy next Full Face NV

CE 0197

 **Manufacturer**  
 Löwenstein Medical  
 Technology GmbH + Co.KG  
 Kronsaalweg 40  
 22525 Hamburg, Germany  
 T: +49 40 54702-0  
 F: +49 40 54702-461  
[www.loewensteinmedical.com](http://www.loewensteinmedical.com)

## Reprocessing Breathing masks


LMT 67506h

**LÖWENSTEIN**  
 medical


**LÖWENSTEIN**  
 medical

Nasal


JULIA

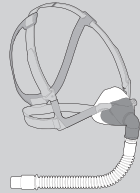
CARA

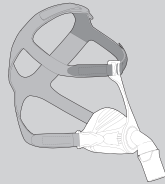JOYCEone

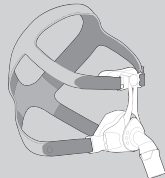
JOYCEeasy

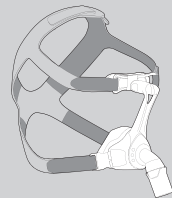
Full Face


LENA

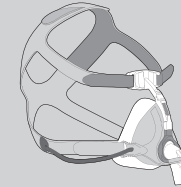
JOYCEone  
Full Face NV

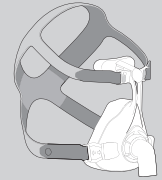
LENA NV

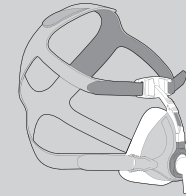
JOYCEeasy  
Full Face

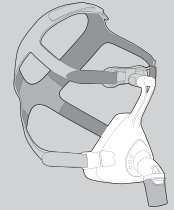
CARA  
Full Face

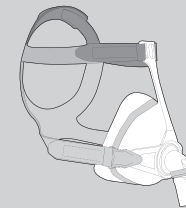
JOYCEeasy next  
Full Face

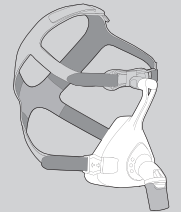JOYCEone  
Full Face

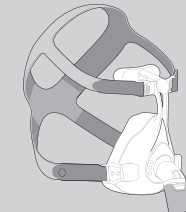JOYCEeasy next  
Full Face NV
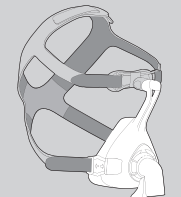
## 1 Introduction

These instructions describe reprocessing in the hospital environment for all Löwenstein Medical Technology mask types and mask parts.

This document is based on standard EN ISO 17664-1 *Processing of health care products. Information to be provided by the medical device manufacturer for the processing of medical devices Part 1: Critical and semi-critical medical devices* (DIN EN ISO 17664-1:2021).

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## 2 Important information

- Follow the information in the respective instructions for use and safety data sheets.
- If the person performing reprocessing (in the hospital) uses a different reprocessing method or exceeds the number of disinfection cycles, the person performing reprocessing accepts responsibility for the resulting effect on the masks.
- National authorities may permit or demand the use of an alternative reprocessing method. In that case, this method must be validated by the person performing the reprocessing.

## 3 Disinfectants

### NOTE

#### Damage to the mask as a result of using an unapproved disinfectant or using a disinfectant incorrectly!

Use of an unapproved disinfectant, incorrect solution strength or disinfection at a higher temperature or for a longer time to take effect than recommended may damage or discolor mask parts and invalidate the warranty for the mask.

- ⇒ Do not use alcohol-based disinfectant.
- ⇒ Use only the cleaning and disinfection methods listed in these instructions.

### 3.1 Disinfectants for chemical disinfection

The disinfectants listed are suitable in terms of their material compatibility and efficacy. Evidence of compatibility with Löwenstein Medical Technology masks has been provided by an independent, accredited test laboratory.

For Germany: The details are based on the list of disinfectants from the Verbund für Angewandte Hygiene e.V. [VAH - German Association for Applied Hygiene].

Disinfectant	Active ingredient basis	Individual active ingredients
Sekusept® Aktiv	Peroxide compounds	Reaction product of peracetic acid which corresponds to sodium percarbonate, non-ionic surfactants, and phosphates

### CAUTION

#### Risk of injury due to unsuitable disinfectants!

Aldehyde-based disinfectants (e.g. gigasept® FF (new)) are not suitable for reprocessing masks of the manufacturer Löwenstein Medical Technology.

⇒ Do not use aldehyde-based disinfectants.

### 3.2 Disinfectants for thermal disinfection

Use according to the manufacturer's instructions

Disinfectant	Type of disinfectant	Individual active ingredients
Neodisher Z	Neutralizing agent	Phosphoric acid, citric acid
Neodisher MediClean forte	Mild alkaline detergent	Non-ionic and anionic surfactants, enzymes

## 4 Reprocessing methods

These reprocessing methods are validated by Löwenstein Medical Technology and permitted for the following masks:

JULIA<sup>1)</sup>, CARA, CARA Full Face, LENA, LENA NV, JOYCEone<sup>2)</sup>, JOYCEone Full Face<sup>2)</sup>, JOYCEone Full Face NV<sup>2)</sup>, JOYCEeasy, JOYCEeasy Full Face, JOYCEeasy next Full Face, JOYCEeasy next Full Face NV

Material of mask parts	Chemical disinfection or thermal disinfection (Vario TD)	Thermal disinfection (autoclave)
Plastic	Max. 30 cycles <sup>4)</sup>	Not permitted
Silicone	Max. 30 cycles <sup>4)</sup>	Max. 30 cycles <sup>4)</sup>
Textile <sup>3)</sup>	Not permitted	Not permitted

<sup>1)</sup> Only subject mask tube to chemical disinfection or replace it.

<sup>2)</sup> It is not necessary to remove the forehead support including the silicone springs for cleaning.

<sup>3)</sup> Replace the headgear on change of patient.

<sup>4)</sup> **Caution!** Exceeding the specified reprocessing cycles can lead to skin irritation.

## 5 Performing reprocessing

You should perform an additional disinfection step after every cleaning step, particularly on change of patient. You can select either a chemical or a thermal disinfection method for disinfection.

### 5.1 Cleaning equipment for manual cleaning and chemical disinfection

For the best possible cleaning result, follow the instructions from the manufacturer of the detergents and disinfectants. Adhere to the quoted concentrations and times to take effect.

Have the following cleaning equipment to hand:

- Solution for washing: 1 ml mild detergent in 1 l water (30 °C - 40 °C). (Mild detergents contain no aggressive constituents and do not irritate the skin. Furthermore, they do not contain perfumes, oils or alcohol.)
- Solution for disinfection: Sekusept® Aktiv, concentration 2%. (A bactericidal, levurocidal, tuberculocidal and mycobactericidal effect is achieved at this concentration and a time to take effect of 15 minutes.)
- Soft cleaning brush
- Soft cleaning brush for mask tube: maximum diameter 10 mm, minimum length 200 mm
- Soft cleaning brush for anti-asphyxia valve: maximum diameter 10 mm, minimum length 100 mm
- Clean tap water (at least drinking water quality)
- Sterile distilled water or fully demineralized water

### 5.2 Cleaning the mask and subjecting it to chemical disinfection

1. Follow infection prevention and health and safety at work regulations.
2. Dismantle mask (see instructions for use). Place all mask parts in a basket or similar to prevent different masks becoming confused.
3. Clean the mask by hand using the solution for washing in accordance with the table below:

Mask part	Action
Mask tube (JULIA only)	Soak for 15 minutes and wash for 3 minutes. Remove the mask tube half-full of solution for washing and take hold of both ends. Move the solution for washing to and fro inside the mask tube for 1 minute. Use the mask tube brush to clean the inside of the mask tube for 3 minutes.
Headgear Frame with frame cover and headbands (JULIA only)	Wash for 15 minutes
Other mask parts	Soak for 15 minutes and wash for 3 minutes. Clean for 3 minutes with a cleaning brush, paying attention to all folds and cavities. Rotate and swivel rotating parts of the mask which cannot be removed (rotating joint) at least 5 times.
Anti-asphyxia valve (full-face masks only)	Clean the openings in the anti-asphyxia valve for 3 minutes using the narrow brush. Carefully lift the inner valve diaphragm and clean it for 3 minutes. Clean both sides of the bearing web.

4. After cleaning: Rinse all parts in clean tap water for 1 minute.
5. Clean all parts in disinfectant solution in accordance with the table (see step 3).
6. After disinfection: Rinse all parts in sterile distilled water or fully demineralized water for 1 minute.
7. Allow all parts to air-dry. If present: Hang the mask tube up to dry.
8. Carry out a visual inspection for cracks and deformations. Replace damaged parts. The silicone mask parts may be discolored or emit a slight odor - this is safe.
9. Remove all parts from the basket or similar and reassemble the mask (see instructions for use).

## 5.3 Subject the mask to thermal disinfection

### Thermal disinfection system

When disinfecting with a certified thermal disinfection system (for Germany: KRINKO/ BfArM recommendations (pages 1252/1254)), adhere to the following time/temperature combination: 90 °C to 94 °C, 5 minutes (range of efficacy in Robert Koch Institute program: A/B).

Observe national and special owner/operator-specific requirements as well as national and international standards and guidelines.

For example:

Miele washer-disinfector running the Vario TD hygiene program with the following water quality, disinfectants, concentrations and program times. (A bactericidal, levurocidal, tuberculocidal, mycobactericidal, sporicidal, fungicidal (clean cond.) and virucidal effect is achieved at this concentration and time to take effect.)

- Dr. Weigert neodisher® Z neutralizing agent at a concentration of 0.1%,
- Dr. Weigert neodisher® MediClean forte universal detergent at a concentration of 0.5%.
- Prewash and main wash: approx. 15 minutes
- Neutralize and rinse: 10 minutes
- Disinfect (total): Disinfect for 25 minutes
- (Tmax): 5 minutes
- Dry: 40 minutes.

Use only demineralized water with a conductivity of ≤ 15µS/cm.

An A0 value of at least 3000 must be achieved.

### Autoclave

Only perform disinfection using an autoclave on change of patient.

Adhere to the following time/temperature combination: 134 °C, 5 minutes, (range of efficacy in a Robert Koch Institute program: A/B).

Observe national and special owner/operator-specific requirements as well as national and international standards and guidelines.

## 6 Packing, storing, and marking

### 6.1 Packing and storing

Packing must stop the mask becoming contaminated with any germ capable of reproducing and pathogenic to humans for the time from after it has been reprocessed until it is used. Use protective packaging to prevent contamination and, if necessary, outer packaging for storage and transport.

Store the mask in its packaging in a dry location away from UV light and dust.

### 6.2 Marking

Mark the reprocessed mask or its packaging as follows:

- Date of reprocessing
- Name of person performing reprocessing

The following markings may also be required:

- Release status (released or not released)
- Address of the person performing reprocessing
- Reprocessing method used

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