
EN-US Instructions for use PL Instrukcja używania HU Használati utasítás CS Návod k použití SK Návod na použitie RO Instrucțiuni de utilizare HR Upute za upotrebu BG Упътване за употреба SL Navodila za uporabo SR Uputstvo za upotrebu LV Lietošanas instrukcija ET Kasutusjuhend LT Naudojimo instrukcija RU Инструкция по пользованию MK Патство за употреба TH คู่มือการใช้งาน MS Arahan penggunaan VI Hướng dẫn sử dụng

WM 65861g 04/2024-EN-US, PL, HU, CS, SK, RO, HR, BG, SL, SR, LV, ET, LT, RU, UK, TH, MS, VI

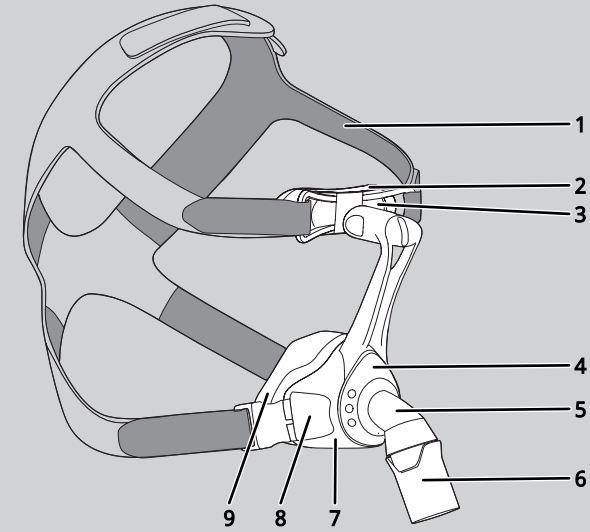CE 0197

 **Manufacturer**  
Löwenstein Medical  
Technology GmbH + Co.KG  
Kronsaalweg 40  
22525 Hamburg, Germany  
T: +49 40 54702-0  
F: +49 40 54702-461  
[www.loewensteinmedical.com](http://www.loewensteinmedical.com)







**JOYCEeasy**  
Nasal Mask

WM 65861g

**LÖWENSTEIN**  
medical

**LÖWENSTEIN**  
medical

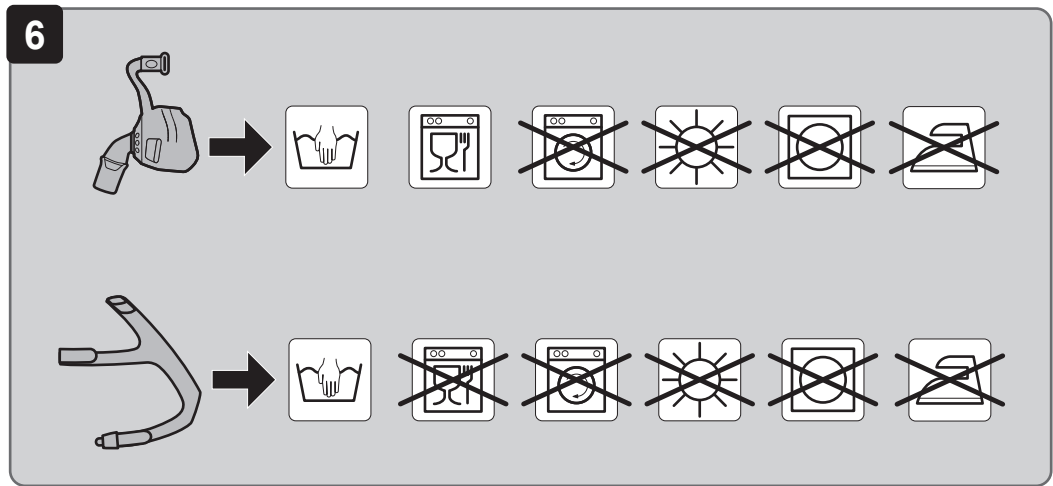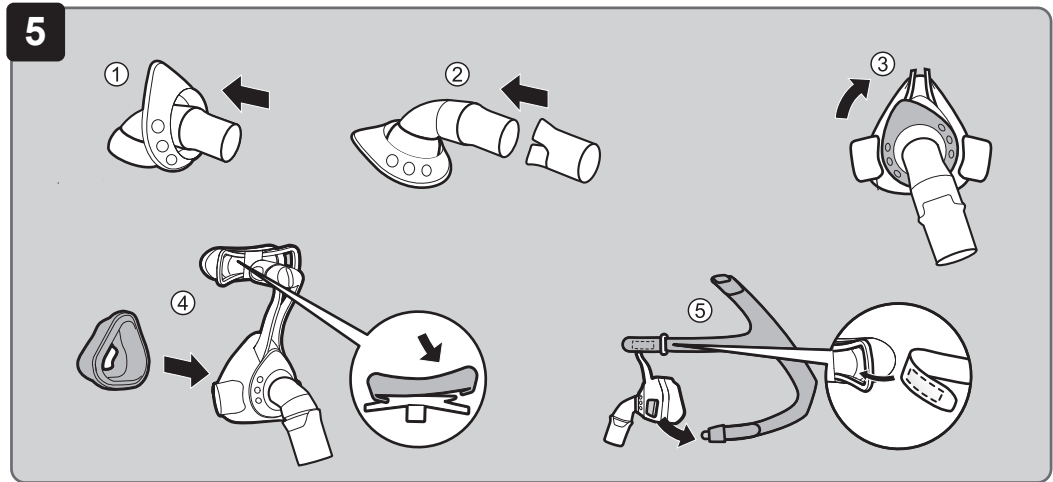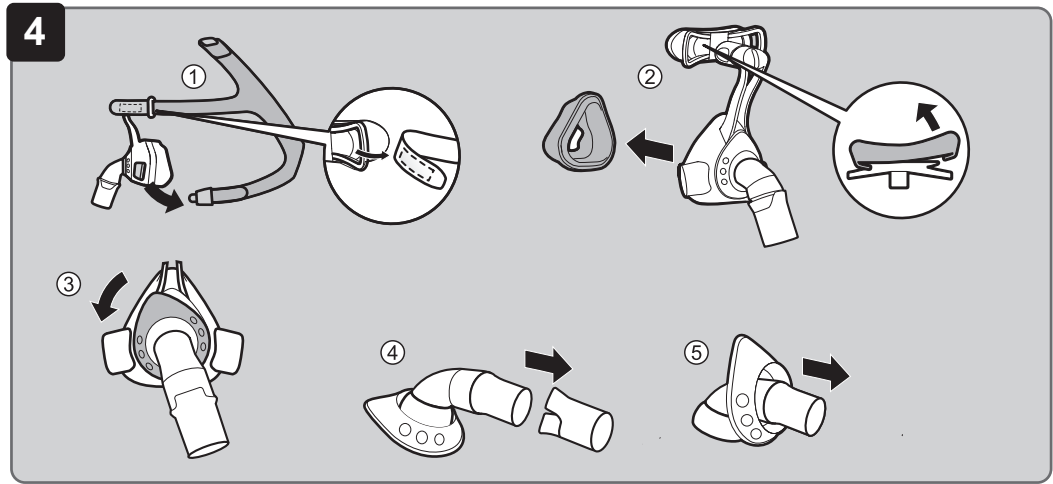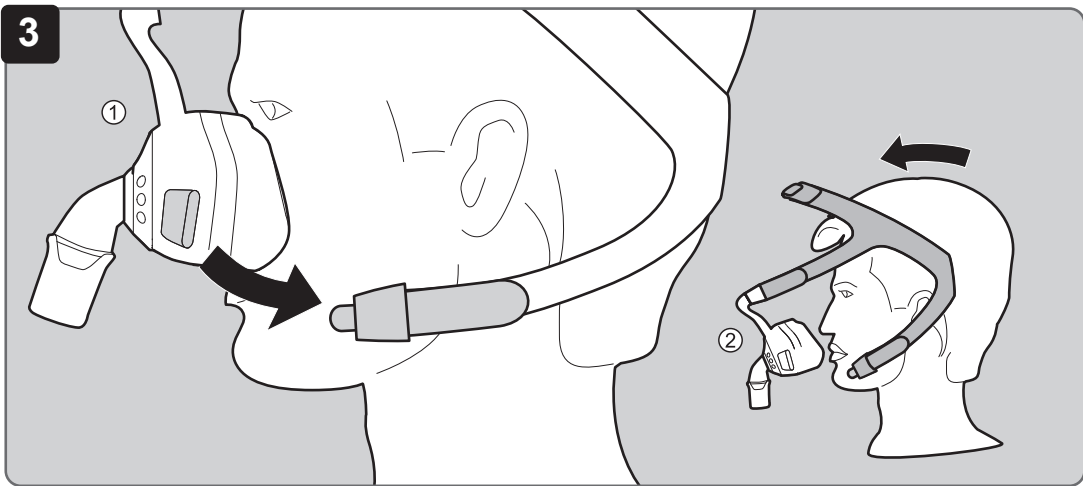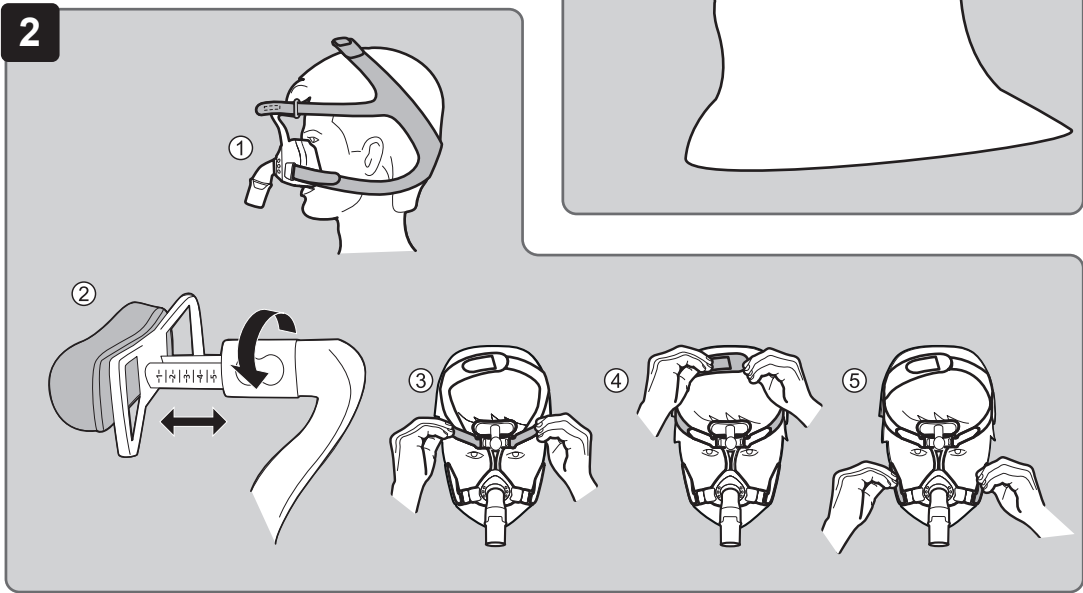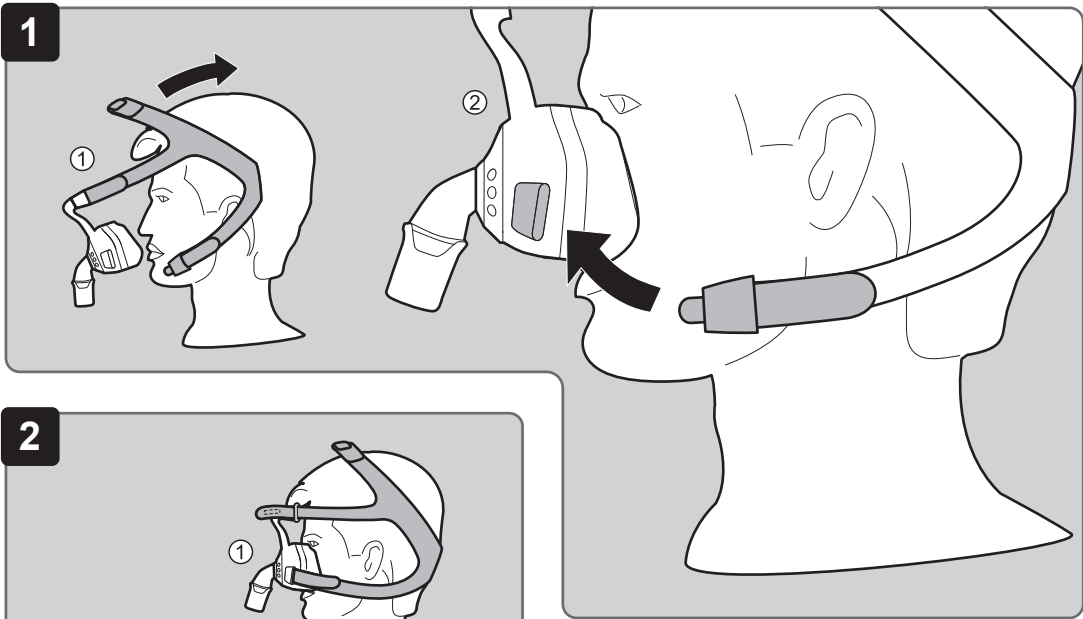with the warranty periods from date of purchase listed below. The warranty conditions are available on the manufacturer's website. We will also send you the warranty conditions on request.

Please bear in mind that any claim under warranty and liability shall be void if neither the accessories recommended in the instructions for use nor genuine spare parts are used.

In the event of a claim under warranty, contact your specialist dealer.

Product	Warranty periods
Masks including accessories	6 months

## 12 Declaration of Conformity

The manufacturer Löwenstein Medical Technology GmbH + Co. KG (Kronsaalweg 40, 22525 Hamburg, Germany) hereby declares that the product complies with the relevant provisions of the Medical Device Regulations (EU) 2017/745. The unabridged text of the Declaration of Conformity can be found on the manufacturer's website.

In the EU: As a user and/or patient, you must report any serious incidents occurring in conjunction with the product to the manufacturer and to the responsible authority.

## 1 Operation

The figures show the following steps for operating the mask:

- 1 Putting on the mask
- 2 Adjusting the mask
- 3 Removing the mask
- 4 Dismantling the mask
- 5 Assembling the mask

**i For blind and partially-sighted users**  
An electronic version of the instructions for use is also available on the manufacturer's website.

## 2 Introduction

### 2.1 Intended use

The JOYCEeasy mask is used for treating sleep apnea and for non-invasive and non-life-sustaining ventilation of patients with respiratory insufficiency. It serves as a connecting element between the patient and the therapy device.

### 2.2 Contraindications

The mask must not be used on patients weighing < 30 kg.

The mask must not be used in the following situations: Immediate intubation required, loss of consciousness, acute vomiting.

The mask may be used in the following situations only with particular caution: Pressure points and acute injuries to the skin of the face; skin allergies involving the face; deformities of the face or nasopharynx; acute pain affecting the face; cough reflex restricted or absent; claustrophobia; acute nausea.

If you are not sure whether one of these situations applies to you, consult your healthcare professional. Observe the contraindications in the instructions for use of your device.

### 2.3 Side effects

The following side effects may occur with use of the mask: Nasal congestion, dry nose, dry mouth in the morning, feeling of pressure in the sinuses, irritated conjunctiva, skin rashes, pressure marks on the face, irritating noises when breathing.

If these side effects occur, contact your healthcare professional.

### 2.4 Clinical benefit

Transfers the therapeutic efficacy of the therapy device to the patient.

## 3 Safety

### Risk of injury due to damaged mask parts or those under strain!

- ⇒ Perform a visual inspection before every use and after every cleaning operation.
- ⇒ Note useful life (see section entitled "Technical Data").
- ⇒ Replace mask parts if necessary.

### Risk of injury due to the use of oxygen!

Oxygen can become deposited in clothing, bedlinen, and hair. Supplying oxygen without a safety device can lead to fire.

- ⇒ Use an oxygen safety valve.
- ⇒ Follow the instructions for use for the oxygen supply system.
- ⇒ Set up oxygen sources at a distance > 1 m from the device.
- ⇒ Do not smoke.
- ⇒ Avoid naked flames.
- ⇒ Ventilate the room well.
- ⇒ Keep mask free from oil and grease.

### Risk of injury due to patient receiving inadequate supply!

- ⇒ Activate low pressure/leakage alarms on the ventilator.
- ⇒ Use the appropriate mask size and check that it is securely in position.
- ⇒ Monitor patients with restricted spontaneous breathing.

### Risk of injury from re-inhalation of CO<sub>2</sub>!

- ⇒ Only use the mask when therapy is in progress.
- ⇒ Only use the mask within the quoted therapy pressure range.
- ⇒ Patients unable to remove the mask themselves must be monitored by a nurse.
- ⇒ Do not close off exhalation systems.

### Risk of injury due to escape of anesthetic gas or atomization of drugs!

- ⇒ Do not use the mask during anesthesia.
- ⇒ Do not use the mask to atomize drugs.

### Risk of injury from inadequate cleaning!

- ⇒ Clean mask parts before using for the first time (see section entitled "Cleaning and reprocessing").
- ⇒ Clean the mask regularly.
- ⇒ When selecting a detergent, consider potential allergies.
- ⇒ On change of patient in a hospital environment: Comply with the document entitled *Information on reprocessing* (see section entitled "Reprocessing").

⇒ For patients with a compromised immune system or particular background of illness, disinfect mask parts daily following consultation with the healthcare professional.

## 4 Product description

A diagram of the individual parts can be found on the title page.

1	Headgear	6	Rotating sleeve
2	Forehead cushion	7	Mask body
3	Forehead support	8	Headgear clip
4	Retaining ring	9	Mask cushion
5	Elbow	10	

### Compatible devices

In some device combinations, the actual pressure does not correspond to the therapy pressure displayed by the device. Have the device adjusted by a healthcare professional so that the actual pressure in the mask corresponds to the therapy pressure. This adjustment should be made with the type of mask that will be used during therapy.

### Exhalation system

The mask has an integrated exhalation system. The exhaled air escapes through a gap.

## 5 Cleaning and reprocessing

### 5.1 Clean mask

1. Wash your hands before starting cleaning.
2. Dismantle mask (see Figure 4).
3. Clean the mask by hand (max. 30 °C, 1 ml mild detergent to 1 l water) in accordance with the following table:

Mask part	Frequency	Action
All mask parts	Daily	Soak and wash for 15 minutes and clean for 3 minutes using a soft cleaning brush.
Headgear	Weekly	Wash for 15 minutes.

**i** All parts can be washed in a dishwasher once a week (max. 70 °C, mild detergent, max. program length 90 minutes, top basket, separate rinse).

4. Rinse all parts with clean water.
5. Allow all parts to air-dry.
6. Carry out a visual inspection for cracks and deformations. Replace damaged parts. Discoloration is harmless.

7. Re-assemble mask (see Figure 5).

### 5.2 Reprocessing (clinical environment)

In the event of a change of patient, follow the instructions in the document entitled *Information on reprocessing*. The document can be found on the manufacturer's website. We will send you the document on request.

### 5.3 Disposal

Dispose of the mask in domestic waste. In the clinical environment: Dispose of the mask in accordance with hospital regulations.

## 6 Troubleshooting

Fault	Cause	Action
Pain due to pressure on the face	Mask too tight.	Loosen headgear.
Draft in the eye	Mask too loose.	Tighten headgear.
	Mask does not fit.	Contact your specialist dealer.
Therapy pressure not reached.	Mask not correctly adjusted.	Re-adjust mask.
	Mask cushion damaged.	Replace mask cushion.
	Patient circuit damaged.	Check circuit and correct fit of circuit.

## 7 Technical data

Dimensions in mm (W x H x D)	
Size S	100 x 100 x 80
Size M	105 x 110 x 80
Size L	110 x 120 x 80
Weight	
Size S	94 g
Size M	97 g
Size L	98 g
Dead space	
Size S	100 ml
Size M	105 ml
Size L	115 ml
Tube connection: Tapered connection to EN ISO 5356-1	Ø 22 mm (male)

Flow resistance at 50 l/min	0.3 hPa
at 100 l/min	1.0 hPa
Service life	5 years
Useful life	Up to 12 months <sup>1</sup>
Therapy pressure	4 hPa - 30 hPa
Quoted two-figure noise emission value according to ISO 4871:	
Sound pressure level	17 dB(A)
Sound power level	25 dB(A)
Uncertainty factor	3 dB(A)
Temperature:	
Operation	+5 °C to +40 °C
Transport and storage	-20 °C to +70 °C
Standards applied	EN ISO 17510: 2020
Product class to MDR (EU) 2017/745	IIa


<sup>1</sup> Useful life depends on cleaning and on the detergents used, on the amount of time worn daily, on therapy pressure, and on individual secretion of sweat. Using the second mask cushion (included as an option) does not extend useful life.

## 8 Materials

No parts of the mask contain latex, PVC (polyvinyl chloride) or DEHP (diethylhexyl phthalate). If you are allergic to any of the materials listed, only use the mask after consulting a healthcare professional.

Headgear	CO (cotton), EL (elastane), P (polyester), PU (polyurethane)
Headgear clip	PA (polyamide), POM (polyoxymethylene)
Forehead support	PA (polyamide)
Forehead cushion	SI (silicone)
Retaining ring (hard component)	PP (polypropylene)
Retaining ring (soft component)	TPE (thermoplastic elastomer)
Mask body	PA (polyamide)
Mask cushion	SI (silicone)
Elbow	PA (polyamide)
Rotating sleeve	PA (polyamide)

## 9 Characteristic pressure/flow curve
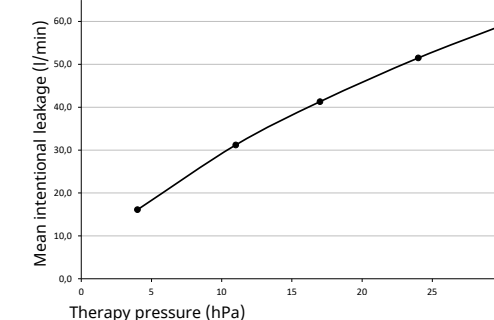
## 10 Markings and symbols

The following markings and symbols may be applied to the device, accessories or packaging.

Symbol	Description
	Manufacturer and, if necessary, date of manufacture
	Unique device identifier (uniform product code for medical devices)
	Order number
	Indicates the product is a medical device
	Follow the instructions for use
	Permitted temperature range for transport and storage
	Use by date
	Keep out of sunlight
	Lot number
	CE symbol (confirms that the product conforms to the applicable European directives/regulations)

## 11 Warranty

Löwenstein Medical Technology gives the purchaser of a new original Löwenstein Medical Technology product and of a spare part fitted by Löwenstein Medical Technology a limited manufacturer warranty in accordance with the warranty conditions applicable to the product in question and in accordance

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